



MIAMI LUMBAR® TLSO. BUILDS OFF THE MIAMI LUMBAR UNIVERSAL BELT TO DELIVER SUPERIOR IMMOBILIZATION

One for all.

Building on the modularity of the Miami Lumbar LSO, the TLSO Kit helps you reduce your inventory, while enabling you to properly and confidently fit nearly every patient you see. The Miami Lumbar universal belt, a selection of panels, and a T-Connector Kit are all you need on-hand.



Easy to fit.

There are only 3 simple steps necessary to fit the Miami Lumbar TLSO. First fit the Miami Lumbar LSO as you normally would. Then set the desired height of the anterior and posterior thoracic extensions. Finally, adjust the twin-angle hinges and cinch down the trimmable shoulder straps.

Successful outcomes.

As your patient progresses through different levels of healing, you can step-down the TLSO to a standard LSO brace, or remove the ATE and add the optional auxiliary shoulder straps to the PTE. Whatever level of immobilization you need to achieve, the Miami Lumbar TLSO accommodates you at every step.



ORDERING INFORMATION

| Part Number | Description |
|-------------|----------------------------|
| LMB-T01 | TLSO system with ATE & PTE |
| LMB-T02 | PTE only |
| LMB-A01 | Air Bladder |
| LMB-A02 | Auxiliary Straps |

ANATOMY OF THE MIAMI LUMBAR™ TLSO



Attachment site with slot secures posterior panels to TLSO belt loop (see below)




TLSO belt loop fits into attachment slot on posterior panels (see above) for increased stability



Optional auxiliary straps for the highest level of immobilization



Adjustable air bladder to compensate for post-surgical changes in lordotic curve

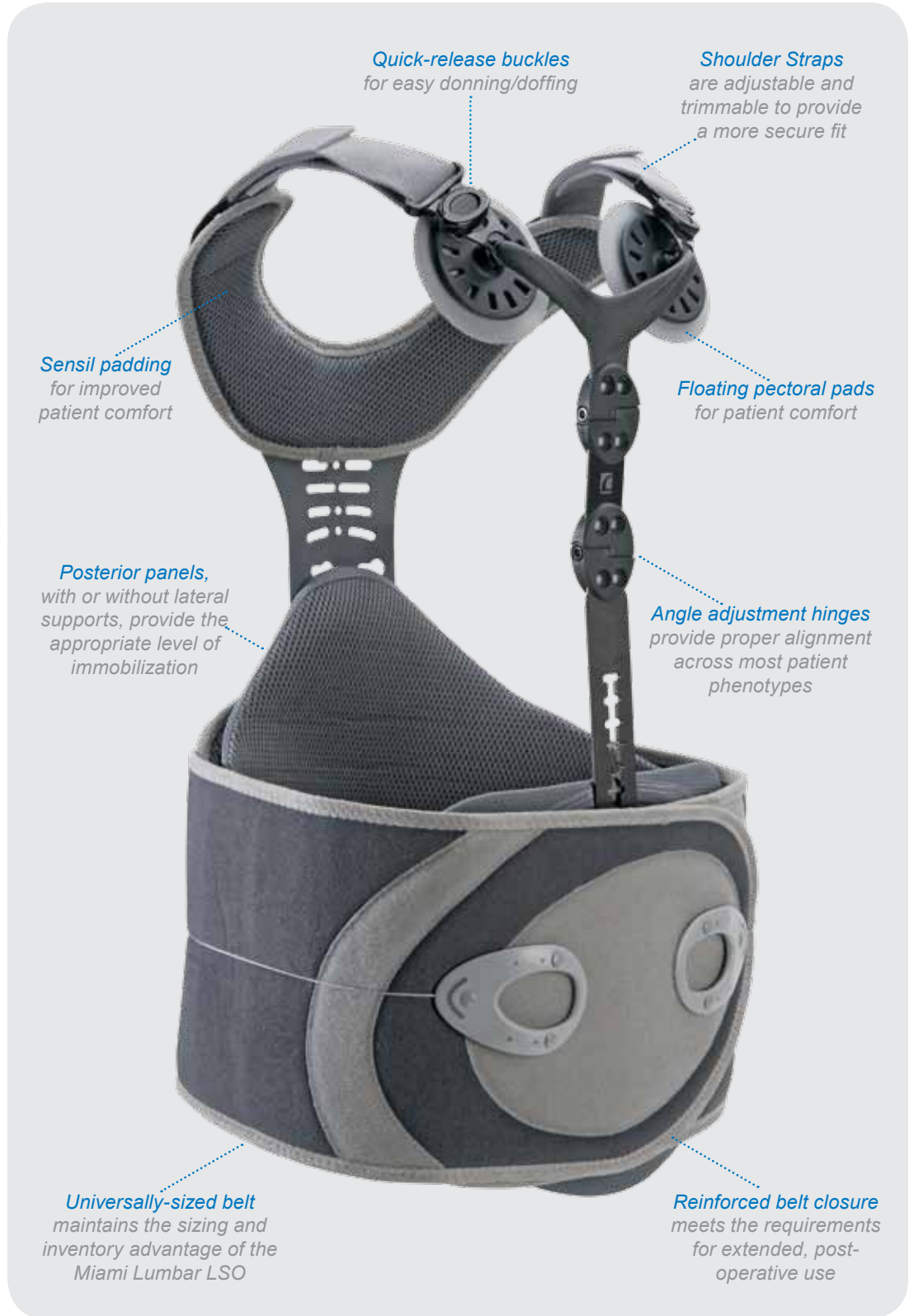

Quick-release buckles for easy donning/doffing

Shoulder Straps are adjustable and trimmable to provide a more secure fit

Sensil padding for improved patient comfort

Floating pectoral pads for patient comfort

Posterior panels, with or without lateral supports, provide the appropriate level of immobilization

Angle adjustment hinges provide proper alignment across most patient phenotypes

Universally-sized belt maintains the sizing and inventory advantage of the Miami Lumbar LSO

Reinforced belt closure meets the requirements for extended, post-operative use

WARRANTY: 6 months

LCODE: L0462

T: 800.767.7776 x 3 • F: 800.869.7776
 estore.spsco.com • www.spsco.com

© 2010 Southern Prosthetic Supply. SPS and Experience Our Commitment are trademarks of SPS. All others are the properties of their respective owners.

SPS
 Experience Our Commitment

OALM-1110R3