


NEW ICEROSS SEAL-IN® X5

PRODUCT FEATURES



IMPACT LEVEL

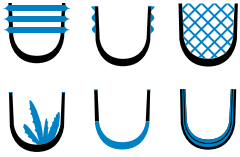

HYPOBARIC SEALING MEMBRANES

Incorporates five serial Hypobaric Membranes which conform to the shape of the internal socket wall, providing an airtight seal.



CREATING HYPOBARIC SUCTION

When the amputee steps into the socket, air is immediately expelled through a distal valve, creating hypobaric suction below the seals.



PREVENTS PISTONING AND ROTATION

Combining the airtight seal with the stabilizing matrix minimizes movement of the limb within the socket.



EASY PRESSURE RELEASE


Depression of the release button on the socket valve allows quick equalization of air pressure for easy removal of the prosthesis



New Iceross Seal-In® X5

New Iceross Seal-In® X5 with Wave

ICEROSS SEAL-IN X5 LINER PROFILE

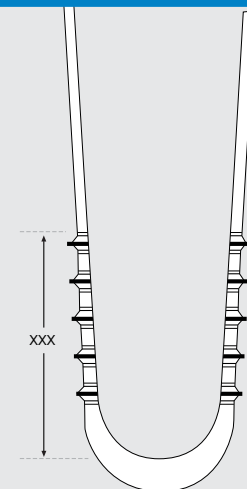

SIZE	XXX(cm)
18	10
20	11
22	11
23.5	11
25	12
26.5	12
28	12
30	12
32	13
34	13
36	13

SIZES

- 11 sizes (18-36) See chart above.

LENGTH

- 300 - 430mm

PROFILE

- 3mm
- Wave 3mm

SEAL TYPE

- Serial (5 Seal Membranes)

MATRIX LENGTH

- 10-13cm, size dependent

PART NUMBERS:

- I-3663XX Iceross Seal-In X5
- I-3664XX Iceross Seal-In X5 with Wave
- XX: 18,20,22,23.5,25,26.5,28,30,32,34,36

6-MONTH WARRANTY