

Silhouette™ VS

- Vertical shock absorption softens impact at heel strike to provide comfort and preserve residual limb health
- Multi-axial function and +/- 15 degrees of inversion/eversion provides excellent ground compliance and stability
- Optimized sole plate with tri-durometer urethane strips transfers load from lateral to medial plane with each step for more natural medial roll-over
- 3/8" or 3/4" heel heights available broaden footwear choices from dress shoes to work boots

Silhouette VS Clearance/Weight Chart

| Size | Clearance 3/8" Heel | Clearance 3/4" Heel | Foot Weight | Shell Weight |
|----------|---------------------|---------------------|-------------|--------------|
| 22-25 cm | 7.2 in. | 7.7 in. | 326g | 127g |
| 26-28 cm | 7.7 in. | 8.2 in. | 391g | 159g |
| 29-31 cm | 8.1 in. | 8.6 in. | 474g | 219g |



+/- 15 degrees of inversion/eversion provides ground compliance and stability for broad range of low to moderate impact K3 level ambulators.

Tri-Durometer strips on carbon fiber sole plate enhances shock absorption and allows medial transfer of ground reaction forces to smooth roll-over.

Indications:

- Low to moderate impact K3 level ambulators
- User Weight Rating: up to 365 lbs

Technical Specifications:

- Sizes: 22-30 (size 31 cm available as special order for additional fee)
- Adapter Options: 34 mm Female Tube Clamp, 34mm Male Tube Clamp, Versa Clamp Kit.
- Footshell: Mid-Profile Shell; High-Profile also available
- Part Number: R14 (3/8" heel); R18 (3/4" heel)


60 Day Trial Period

Warranty: Foot-36 mos.; Footshell - 6mos.

Mfg's Suggested LCode^{††}: L5987, L5986

New Silhouette & Silhouette LP

- Anterior pylon shift of 1 cm (on average) improves heel to toe dynamics.
- NEW 3/4" heel option (Silhouette, RS8 and Silhouette LP, RS9) available.

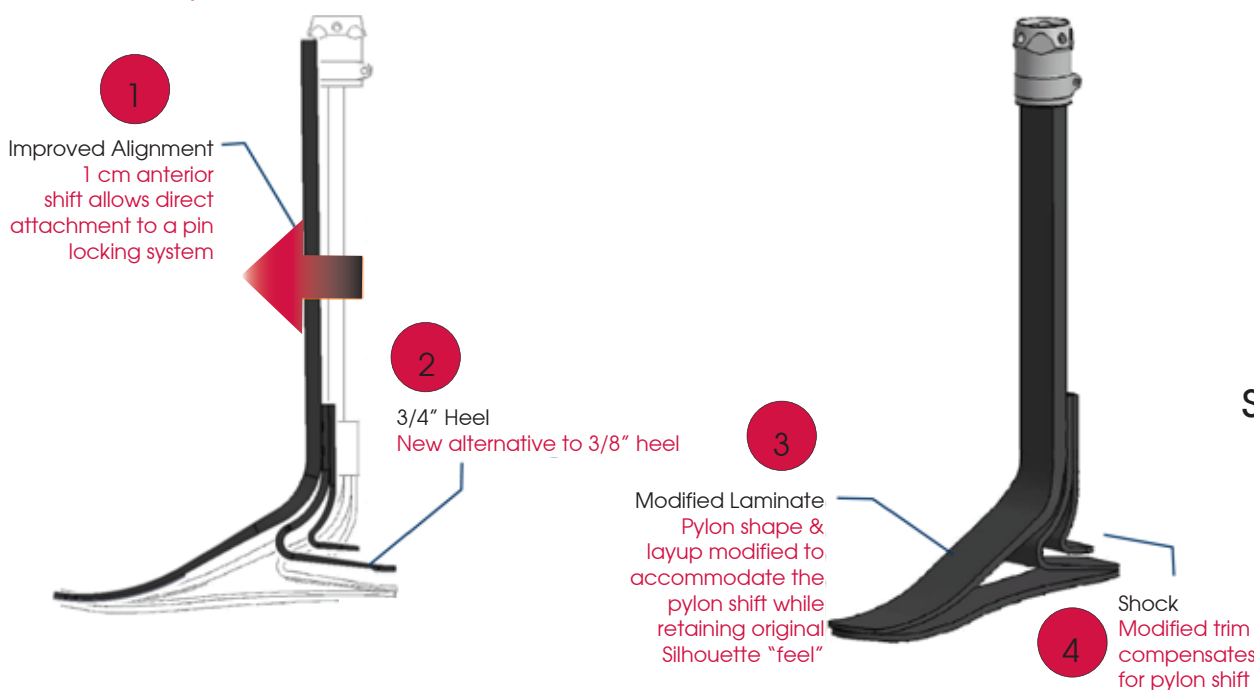

Silhouette (RS4 or RS8)



Silhouette LP (RS5 or RS9)



KEY BENEFITS OF NEW ENHANCEMENTS:

- Direct attachment to a pin locking system simplifies alignment process and expands range of applicable patients
- Optimized heel to toe ratio allows users to experience balanced step and normal stride length
- Easier to align due to anterior pylon shift removes need for a slide adapter under the socket, reducing weight and increasing user comfort & satisfaction
- 3/4" heel height option adds choice and broader range of applicable footwear

NEW Silhouette & Silhouette LP Clearance/Weight Chart

| | Size | Clearance 3/8" Heel | Clearance 3/4" Heel | Foot Weight** | Shell Weight |
|------------------------------|----------|---------------------|---------------------|---------------|--------------|
| NEW & Improved Silhouette | 22-25 cm | 7.2 in. | 7.7 in. | 506g | 127g |
| | 26-28 cm | 7.7 in. | 8.2 in. | 536g | 159g |
| | 29-31 cm | 8.1 in. | 8.6 in. | 566g | 219g |
| NEW & Improved Silhouette LP | 22-25 cm | 6.8 in. | 7.3 in. | 357g | 127g |
| | 26-28 cm | 7.3 in. | 7.8 in. | 387g | 159g |
| | 29-31 cm | 7.7 in. | 8.2 in. | 417g | 219g |

** Foot Weight in chart left includes adapter kit

60 Day Trial Period

Warranty: Foot-36 mos.; Footshell - 6mos.

Mfg's Suggested LCode^{††}: L5980, L5986

^{††}All coding opinions expressed herein are based upon information provided by the manufacturer. Regardless of the source of the information provided, the final responsibility for correct coding within all established laws, rules, standards and practices is the sole responsibility of the practitioner submitting the claim.