

NEW & IMPROVED TOTAL KNEE®

BY ÖSSUR®

The redesigned Total Knee incorporates a full range of improvements:

- Increased durability
- Improved aesthetics
- Easier cosmetic finishing

The Total Knee offers the following advantages:

- Mid-swing shortening provides increased toe clearance to prevent hip-hiking and reduce chances of tripping
- Extension assist and extension promotor ensure full knee extension
- Geometric lock creates stability throughout stance phase
- Adjustable stance flexion bumpers mimic knee flexion during early stance, reducing shock and stress on the residual limb
- Rugged and reliable, yet lightweight

Total Knee 1900 is compatible with the Flex-Foot® Assure.





Total Knee 2000 is compatible with the Flex-Foot® Vari-Flex® with EVO™ or the Talux (far right).



Experience Our Commitment

T: 800.767.7776 x 3 • F: 800.869.7776
estore.spsco.com • www.spsco.com

TECHNICAL SPECIFICATIONS

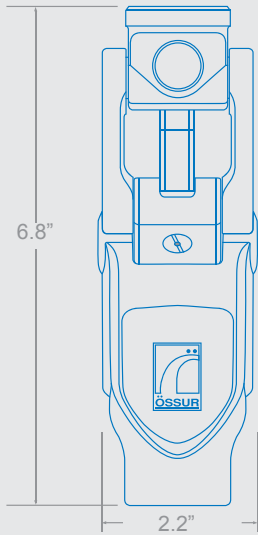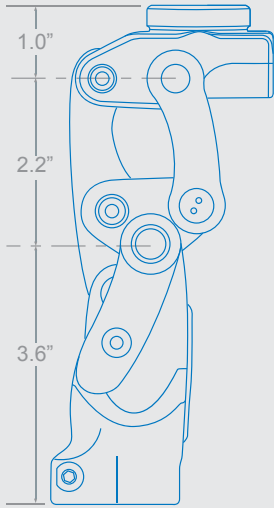

NOTE: Dimensions are identical for both models

TOTAL KNEE® 2000

INDICATIONS:

- Multi speed ambulators who are able, or have the potential to walk at different speeds
- User weight limit of 220 lbs
- Moderate to high impact levels

FEATURES & BENEFITS:

- Polycentric knee with geometric locking system
- Swing control adjustments through 3-valve hydraulic system
- Adjustable stance flexion
- Adjustable extension promoter

TOTAL KNEE® 1900

INDICATIONS:

- Single speed ambulators engaged in low impact activities and constant walking speed
- User weight limit of 220 lbs
- Low to moderate impact levels

FEATURES & BENEFITS:

- Polycentric knee with geometric locking system
- Controlled and steady swing phase through use of adjustable friction
- Adjustable stance flexion
- Adjustable extension promoter

TECHNICAL SPECIFICATIONS:

Identical for 2000 and 1900 except as noted:

- WEIGHT LIMIT: 220 lbs
- PRODUCT WEIGHT: 1.452 lbs (2000); 1.416 lbs (1900)
- STANDARD BUILD HEIGHT: 6 7/8"
- KNEE FLEXION AVAILABLE: 160°
- WARRANTY: 24 mos limited warranty
- TRIAL PERIOD: 30- DAY no obligation

PART NUMBERS:

2000	Total Knee 2000
1900	Total Knee 1900
2070S	Optional Extension Assist

LCODE: xxxx, xxxx