

# Compare Products

800-767-7776 x3

[Contact Your Sales Manager](#)

Compare two or more products and their technical specifications.

Attributes you can compare:

- Pricing
- Description
- Patient weight limits
- Supplier
- L-Codes









# Compare Products

*How it Works*

800-767-7776 x3

[Contact Your Sales Manager](#)

1. Click the  symbol underneath a product you want to compare

	
<div><div>SELECT OPTIONS</div><div></div></div>	<div><div>SELECT OPTIONS</div><div></div></div>
<b>RUN-KIT</b>	<b>ROM-KIT</b>
RUNWAY & RUNWAY HX	KINTERRA FOOT/ANKLE SYSTEM
<i>FREEDOM INNOVATIONS</i>	<i>FREEDOM INNOVATIONS</i>


# Compare Products

## *How it Works*

800-767-7776 x3

[Contact Your Sales Manager](#)

2. The product will add to queue.
3. Add one or more items to compare by repeating step 1

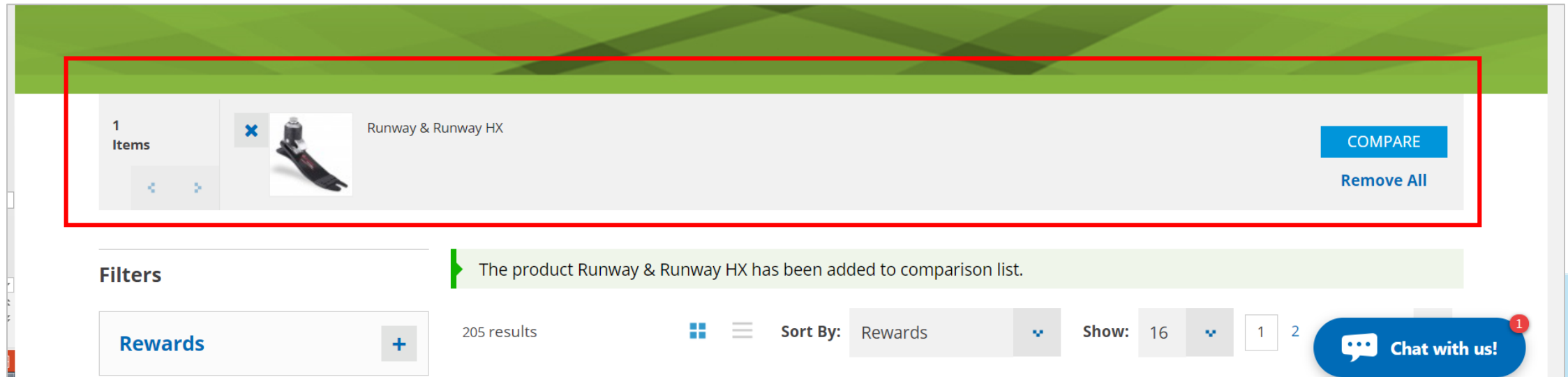


# Compare Products

*How it Works*

800-767-7776 x3

[Contact Your Sales Manager](#)

4. Click **COMPARE** once you've added all the products you'd like to compare

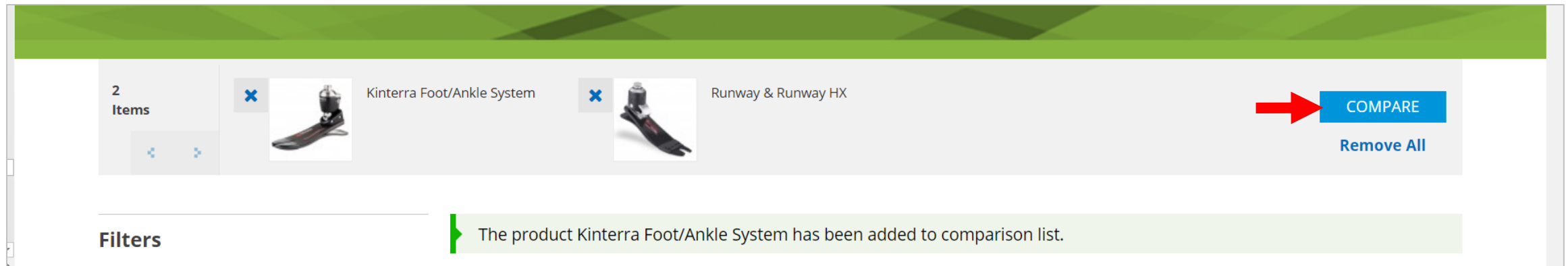

**Pro Tip:** Select no more than 4 products for the best view of all specifications.



# Compare Products

How it Works

800-767-7776 x3


[Contact Your Sales Manager](#)

Print




## 5. Compare specifications



### Compare Products




**RUN-KIT**  
RUNWAY & RUNWAY HX  
FREEDOM INNOVATIONS



**ROM-KIT**  
KINTERRA FOOT/ANKLE SYSTEM  
FREEDOM INNOVATIONS

PRICE		
DESCRIPTION	<p><b>Heel Height Adjustments- On the Go!</b></p> <p>Amputees greatly value the versatility offered by the Runway adjustable heel height prosthetic foot. Its anatomical gliding ankle maintains appropriate foot alignment to provide knee stability and consistent performance at all heel heights. Users can adjust from a flat heel for sandals or walking barefoot at home, up to a 2 inch heel for dress shoes and even cowboy boots. The Runway is the only heel height adjustable product manufactured using carbon fiber, which stores and returns energy so that users can walk farther, faster, longer. Unlike competitive products, the Runway incorporates carbon fiber layers that are tailored to an individual's body weight and impact level, providing a customized product for superior comfort. Rated to 255lbs., the Runway offers the highest user weight rating in the heel height adjustable product category, making it accessible to more amputees.</p> <p>Originally launched with a push button adjustment mechanism with ten height settings, the new Runway HX allows users to make precise adjustments within the 0-2" range. A port in the posterior (rear) aspect of the foot shell provides access to adjust the ankle by simply turning a 4mm hex key to reach the desired position. In addition to the ability to make easy and precise adjustments, the Runway HX also offers tighter tolerances for improved durability.</p> <p>Both the push button Runway and the Runway HX are available with an optional sandal toe foot shell, which has been redesigned to appear more natural and is capable of pleasing even the most discerning users.</p> <ul style="list-style-type: none"><li>As the only user-adjustable heel height foot designed for active K3 amputees, the device is customized to user weight and activity to provide superior comfort and energy return</li><li>Adjust heel height from 0-2" with option of push button or 4mm hex key adjustable design</li><li>Anatomical "gliding" ankle design maintains appropriate foot alignment, regardless of heel height</li><li>Runway HX offers adjustments at any place within the 0-2" range, compared to 10 preset adjustments with the push button design</li><li>New cosmetically improved sandal toe foot shell for improved user acceptance</li></ul>	<p><b>Promotions</b></p> <p><b>Freedom: FREE Foot Shell Promotion</b></p> <p><b>through December, 2017</b></p> <p>Purchase a Maverick or Kinterra and get a FREE foot shell.</p> <div></div> <p><b>Freedom Innovations: Ideal Combo</b></p> <p><b>Extended through December 31, 2017</b></p> <p>For a limited time only, purchase a Plie 3 and get 50% off the Kinterra foot/ankle system (list price)</p> <p>OR</p>

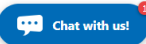

**Pro Tip:** You can print this page.

